

AcceleDent®

OPTIMA™

PRACTICE GUIDE



1st & Only Orthodontic Device that Directly Connects Patients and Practices with Monitoring, Direct Messaging and Rewards

AcceleDent Optima offers a state-of-the-art solution for encouraging patient engagement and assessing compliance when patients leave your office. With the AcceleDent App, you and your practice will have direct and immediate knowledge of compliance, and will be able to message non-compliant patients right away to encourage use. This proactive motivation, when combined with AcceleDent's SoftPulse Technology®, helps you to achieve the clinical outcomes you expect from your treatment plan.

Partner with your Patient

AcceleDent Optima's cloud connection capabilities connect patients to your practice to optimize the treatment experience:



- At-a-glance usage data for your AcceleDent Optima patients, allowing for proactive case management
- Personal information is secured with HIPAA compliant encryption



- Ability to direct message patients to encourage compliance and maximize efficiency of treatment
- Patients can directly communicate back to the practice to feel more "connected" with their treatment team between appointments. They don't have to call or wait, but can send messages through the app.

AcceleDent Optima Doctor Experience Overview



The AcceleDent Practice Account can be downloaded as an app or accessed at MyAcceleDent.com. The user interface is very similar to the patient-facing app, making it easy to navigate and learn.

The Practice Account makes it convenient for the doctor and staff to all view the same account as a team and see all active AcceleDent patients in one place. Key features include:

- Patients are searchable and sortable by name and compliance ranking
- The Ranking shows each patient's overall compliance



- Your practice can focus messaging or encouragement to the patients that need it the most
- Click on a patient's name to be taken to their individual dashboard to see their usage over the past 30 days and details including: days of correct use, overuse and missed sessions in the last 30 days and the time of day the patient used their AcceleDent

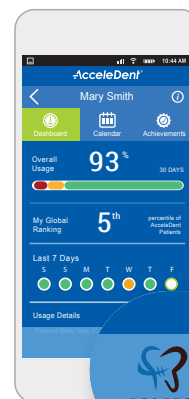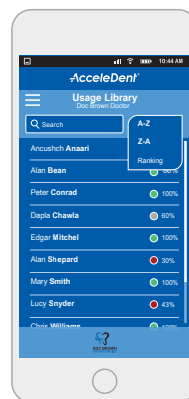

A practice logo may also be added to your practice account. This logo will also show up on the bottom of your patient's app.



Designed to Benefit Your Practice and Patients

AcceleDent Optima combines high-tech performance features and a comprehensive app to enhance your patient's accelerated orthodontic experience, empowering them to become a full participant in their treatment.



Features designed to fit your patient's lifestyle

- Small & Lightweight
- Waterproof Activator
- Wireless Charging
- Bluetooth® Enabled

Patented SoftPulse Technology

- Demonstrated to speed up bone remodeling. Pulsating forces, like those AcceleDent produces, are clinically shown to stimulate cellular activity in orthodontic treatment.
- 15 peer-reviewed articles, including 3 randomized control trials.

AcceleDent App promotes patient cooperation

- Users share progress within their **Community** and treatment team
- **Motivation** tools encourage usage
- **Dashboard** reports track and monitor compliance

The AcceleDent App Patient Experience Overview

The AcceleDent App utilizes well known behavioral science techniques to effectively motivate positive behavior and transform the patient's accelerated orthodontic treatment.

The patient uses the app to check their dashboard for progress toward achievements, share their treatment progress with family and friends, connect with you and your practice for important information and set up daily treatment reminders.

Not only can patients use the AcceleDent App, but your practice and the patient's family and friends can also use the app as Observers to support their accelerated orthodontic treatment journey.

Key AcceleDent App Features

On the **My Usage** view, the patient will be able to view session history in a dashboard view or as a calendar.

On the **Dashboard** view, the patient will be able to see their overall usage percentage, which is their own achievement progress towards the goal of daily 20 minute sessions, their global ranking compared to all the other AcceleDent patients, and more usage details.



On the **Calendar** view, the patient will achieve a green circle for days that they have used their AcceleDent Optima for the recommended 20 minutes, orange for days with more than 20 minutes, and red for days with less than 20 minutes of use. The patient will also be able to click through on each day to see what time of day their treatment session occurred.



On the **Achievements** Screen, the patient will see awards they have earned and awards that are still waiting for them to unlock.



Talk to your OrthoAccel representative about AcceleDent Optima today!

FOR MORE INFORMATION | INQUIRIES IN NORTH AMERICA Call 1-866-866-4919 or visit acceleident.com
INQUIRIES OUTSIDE NORTH AMERICA Call +49-800-724-4532 or visit acceleident.eu



Leader in Accelerated Orthodontics®

1 Data on file. AcceleDent® Optima™ is an orthodontic accessory intended for use during orthodontic treatment. Refer to the Directions for Use with AcceleDent Optima or acceleident.com for warnings and precautions prior to using this product. Federal (U.S.) law restricts this device to sell by or on the order of a licensed doctor. Rx use only for a single patient. Bluetooth® is a registered trademark owned by Bluetooth SIG, Inc. OrthoAccel®, AcceleDent®, and SoftPulse Technology® are registered trademarks and Optima™ is a trademark of OrthoAccel Technologies, Inc.