
# AcceleDent®

## OPTIMA™
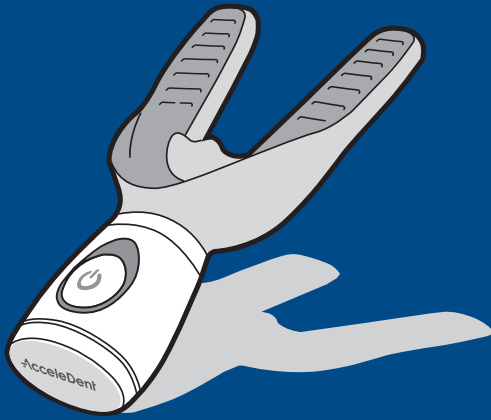
## QUICK REFERENCE GUIDE

REF 405-1302-001 Rev. A


### INTENDED USE/Indications for Use

AcceleDent® Optima™ is an orthodontic accessory intended for use during orthodontic treatment. It is used in conjunction with orthodontic appliances such as braces and helps facilitate minor anterior tooth movement. Its effect on overall treatment time has not been established. This is due to the many individual factors that may impact orthodontic treatment. You are encouraged to speak with your orthodontist to determine the impact this device may have on your specific treatment plan.

### Contraindications for Use


- Use of osteoporosis drugs
- Poor oral hygiene
- Periodontal disease that is not under full control for at least 3-4 months prior to the start of the treatment

# FEATURES AND COMPONENTS

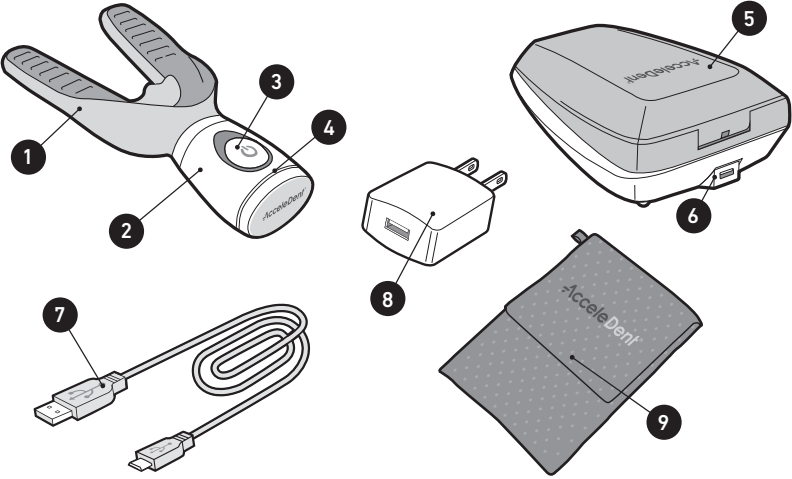
- 1 MOUTHPIECE
- 2 ACTIVATOR
- 3 ON/OFF BUTTON
- 4 STATUS INDICATOR
- 5 CHARGING CASE
- 6 USB CHARGING PORT
- 7 USB CHARGING CABLE
- 8 POWER ADAPTER
- 9 TRAVEL POUCH

# ACCELEDENT USE

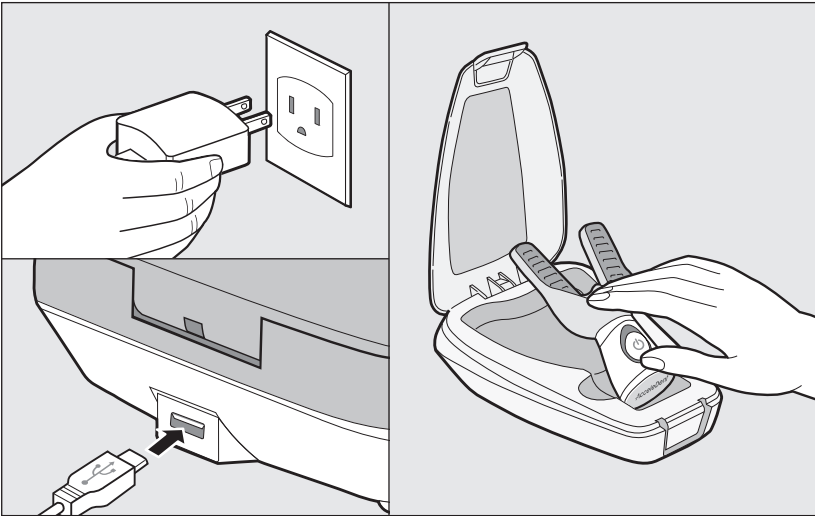
Power the Charging Case by connecting the USB Charging Cable to an electrical source. Prior to first use, place the Activator in the Charging Case to charge for at least 8 hours. The status indicator on the Activator will be solid orange while charging and turn off when fully-charged.

## ACCELEDENT USE (Continued)

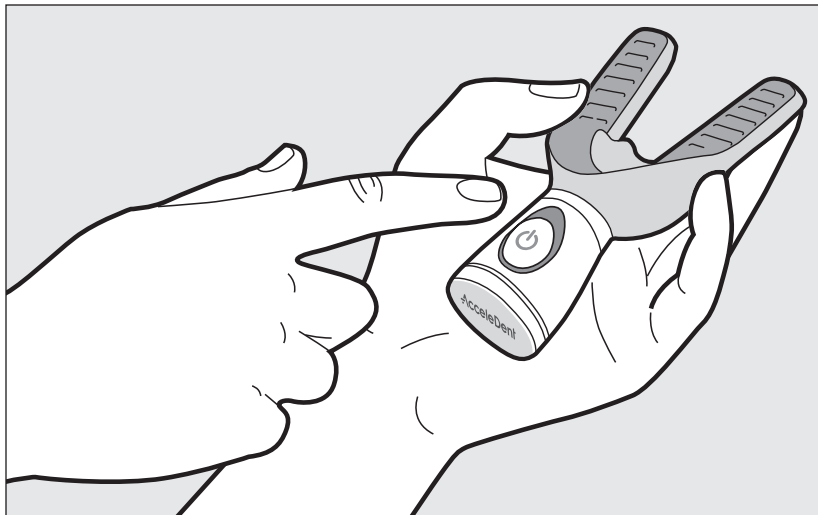
To turn the device on, press the on/off button. The Activator will vibrate and the status indicator will be solid blue.

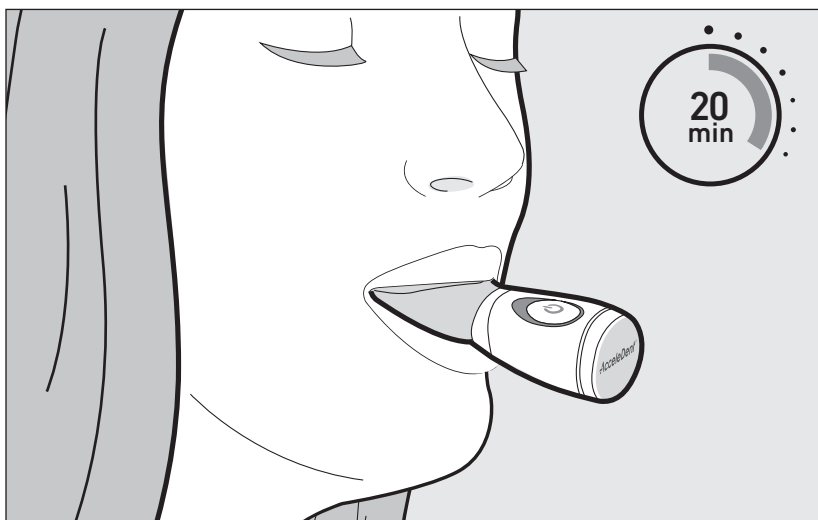Place the mouthpiece in your mouth and bite down gently. After 20 minutes of treatment, the Activator and status indicator light will automatically turn off. To pause your treatment session up to 5 minutes, press the on/off button once. During the pause, the status indicator will pulse blue slowly. The status indicator will blink orange when the battery is low.



## ACCELEDENT USE (Continued)

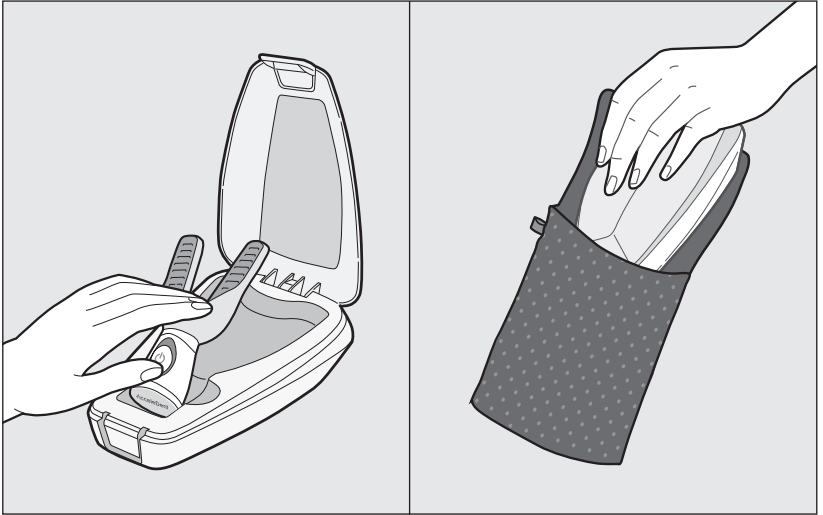Recharge the Activator after each treatment by returning to the Charging Case. A Travel Pouch has been provided to store the Charging Case and accessories when not in use.

## CARE

To clean the Charging Case, unplug from power source and wipe inside and out using a moist cloth and dry with a soft cloth before reconnecting to power source (**DO NOT** rinse or submerge Charging Case).

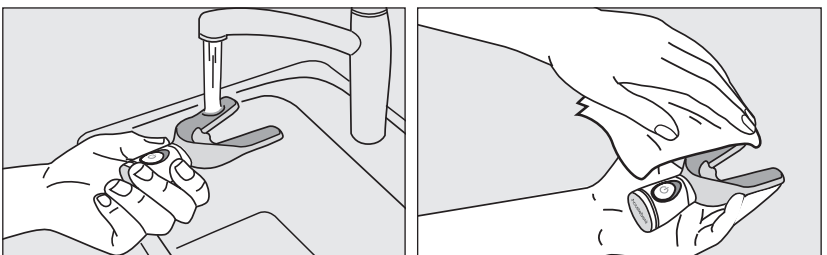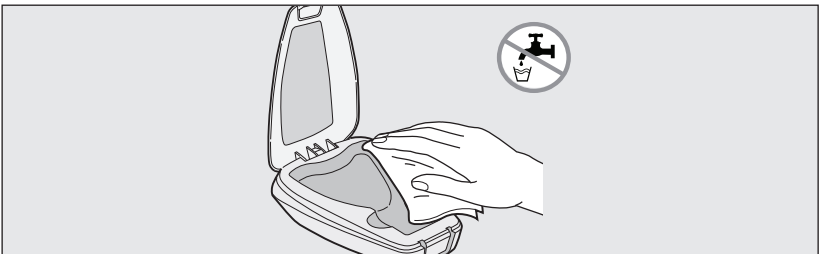Clean the Mouthpiece with warm running water and wipe with a dry cloth prior to storage or charging.

## SAFETY, WARNINGS, PRECAUTIONS

- Use this product only for its intended use as described in the Quick Reference Guide and Directions for Use. Please refer to detailed instructions found on [acceleddent.com](http://acceleddent.com).
- This is a single patient device (one patient to one device).
- This is a prescription device which must be prescribed by an orthodontist or dentist.
- Contact your orthodontist or dentist if any discomfort, or pain occurs during and after use of the product.
- Warning: Do not plug the Power Adapter into an outlet with a voltage other than what is specified on the Power Adapter.
- Warning: Never force the Power Adapter into an outlet.
- Warning: Never use damaged Activator, Charging Case, USB Charging Cable, or Power Adapter.
- Warning: Do not place or store the device or accessories near any heated surfaces.
- Warning: Do not place or store the Power Adapter and Charging Case where it will sit in a pool of water or where it can fall or be pulled into a bathtub, sink, or toilet. Do not reach for a service if this occurs. Charging Case that has fallen into water. Unplug immediately and call customer service if this occurs.
- Do not sterilize the product.
- Do not clean components or accessories in the dishwasher.
- Activator is waterproof, but it is recommended that it is not submerged in water or cleaning solutions for an extended period of time.
- Do not rinse or submerge the entire Charging Case in water.
- Do not attempt to service or perform maintenance on the device. No modification of the equipment is allowed.
- Note: If the LED is alternating orange and blue, the Activator has reached end of life. If you are in the middle of a treatment plan, please contact your orthodontist.

## INDICATORS FOR ACCELEDENT DEVICE

Solid orange	Charging battery
Blinking orange	Low battery
Solid blue	Active treatment
Solid blue and pulsing vibration	Mid-point treatment
Slow pulsing blue	Paused treatment
Fading blue to none	End of treatment
Alternating orange and blue blinking	End of device life

### TO INSTALL THE ACCELEDENT FREE COMPANION APP, PLEASE FOLLOW THE STEPS BELOW:

1. Android™ Device: Search in Android Market or Google Play Store for the AcceleDent App.  
Apple® Device: Search in Apple App Store for the *AcceleDent* App.
2. Follow specific steps recommended by your device manufacture for downloading app.
3. Open the app and begin registration by clicking on SIGN UP.
4. Follow prompts within app for set up and linking the Activator.

\* OrthoAccel is an independent manufacture and has not been authorized, sponsored, or otherwise approved by Apple Inc. or Google Inc. Apple is a registered trademark of Apple Inc. in the U.S. and other countries. Google and Android are trademarks of Google Inc. in the U.S. and other countries.

## CUSTOMER SERVICE

If you have any questions or comments after reading the enclosed materials, please contact OrthoAccel Customer Service:

North America: **1-866-866-4919**


Email: **customerservice@orthoaccel.com**

Online: **www.acceledent.com**

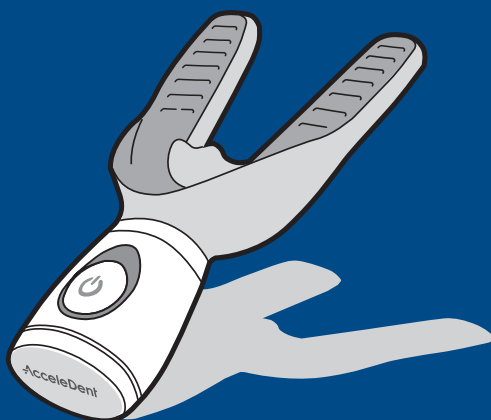Please refer to the AcceleDent Directions for Use for IMPORTANT SAFEGUARDS and detailed instructions found on [acceledent.com](http://acceledent.com).

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